# **REF EASY SHOWER**

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***BRUGSANVISNING***

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Product code: EASY SHOWER

### 1. Introduction

Before using the product check whether the equipment has no transport damage. If any damage is identified, please contact the retail outlet. Then it is important to refer to the content of this manual that provides relevant information for the user. If you have any questions concerning the use of the equipment, please contact us. Our personnel will provide you with the required information on an on-going basis. We would like to point out at the same time that if you have any questions not related to the product use, please contact your doctor, nurse or physical therapist.

### 2. Design

A portable, inflatable basin made of PVC for hair washing with a water drain outlet is soft and comfortable for use. Dual air chambers constitute a protection against spillage of water and suds.


**The kit contains:** basin, water reservoir, shower, drain line, pump, reservoir bracket

### 3. Application


The hair washing kit allows the patient's head to be washed in the prone position. It is particularly recommended for indoor use during care for the bedridden and elderly people with difficulties in tilting their head and bending.

### 4. Method of use


1. First inflate the basin using the pump and close the water drain valve at the bottom of the basin. The basin can be filled with warm water with a temperature below 50°C.

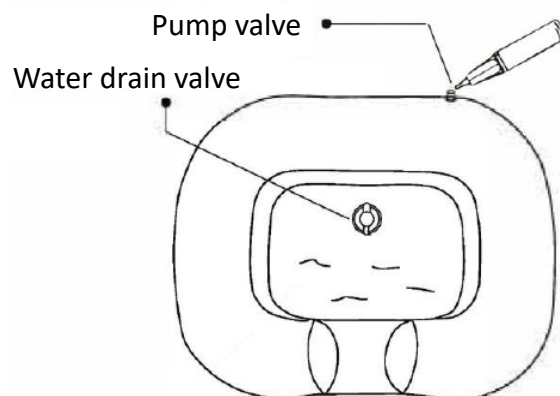
2. Connect the water reservoir to the shower.

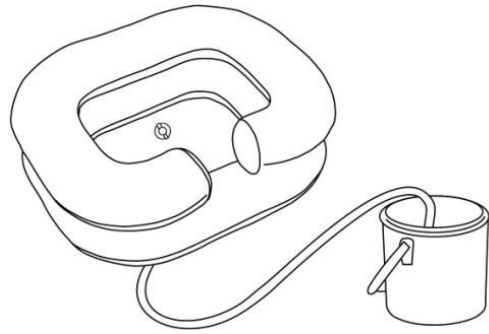3. Fill the reservoir with warm water with a temperature not exceeding 50°C. The reservoir can be hung on a stand or on the wall using the reservoir bracket.

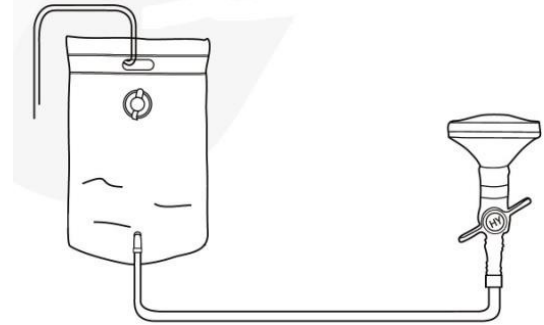
4. Prepare a bucket and place the drain line inside it. To ensure proper water drainage, check that the line is not kinked.

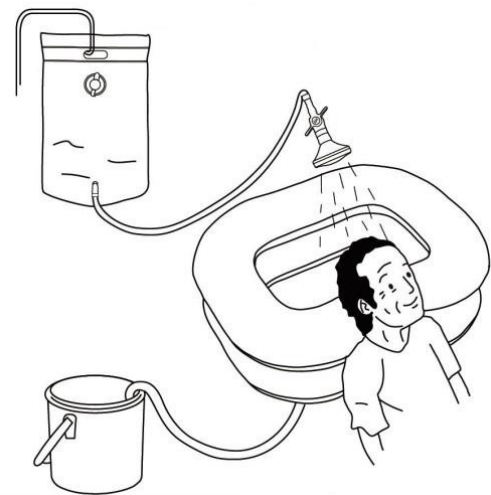5. Place the patient's head in the neck recess and unscrew the shower valve to adjust water flow.

6. After washing the patient's head open the drain valve, so that the water can flow to the bucket through the drain line.
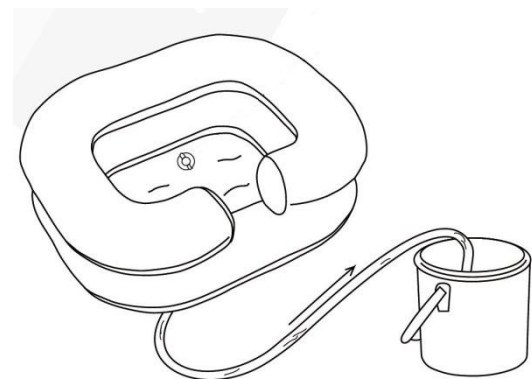
## 5. Technical data

Basin dimensions	65 x 50 x 25 cm
Weight	0,8 kg
Basin volume	10 liters

Water reservoir volume	11 liters
------------------------	-----------

## 6. Notes

**Timago International Group does not bear responsibility of improper use of the basin, failure to comply safety regulations and misuse.**

## 7. Storage and transport conditions

Device should not be exposed to direct sunlight, water or moisture.

## 8. Warranty information

All products distributed by our company are under a warranty, the terms of which are described in the warranty card received at the time of purchase. Please, keep in mind that for warranty purposes, the proof of purchase (receipt or invoice) should be retained.

## 9. Labels

	Reference number		Manufacturer
	LOT number		UDI code
	Serial number		Medical device
	Note		Please read the instruction
	The manufacturer checked for compliance with essential medical devices requirements.		

Our complete line of products is available from our distributors. Learn more at: [www.timago.com](http://www.timago.com). Thank you for choosing Timago!

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